

PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%


Three samples of known concentration were tested twenty times on one plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

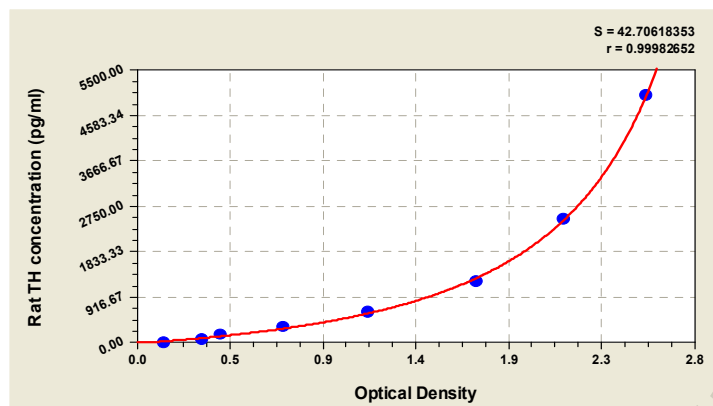

pg/ml	OD1	OD2	Average	Corrected
5000	2.563	2.501	2.532	2.390
2500	2.113	2.131	2.122	1.980
1250	1.683	1.699	1.691	1.549
625	1.142	1.168	1.155	1.013
312	0.712	0.754	0.733	0.591
156	0.418	0.431	0.425	0.283
78	0.322	0.341	0.332	0.190
0	0.141	0.143	0.142	

LOD

19.5 pg/ml.

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of rat TH in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:100	Average %	90
	Range %	85-96
1:200	Average %	96
	Range %	91-105
1:400	Average %	98
	Range %	92-104
1:800	Average %	94
	Range %	86-98

RECOVERY

The recovery of rat TH spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	90	85-94