

PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%


Three samples of known concentration were tested twenty times on one plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

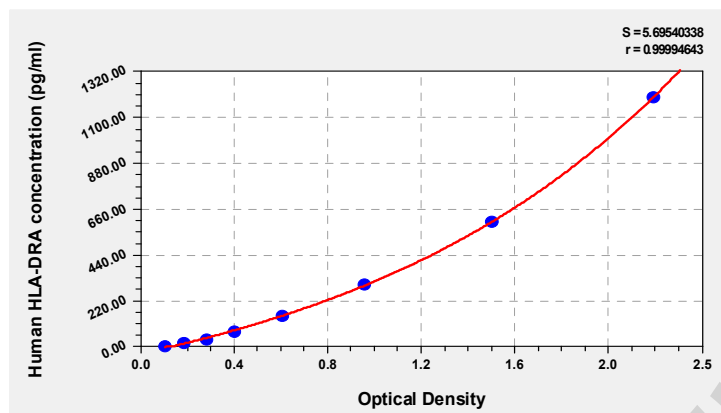

pg/ml	OD1	OD2	Average	Corrected
1200	2.247	2.235	2.241	2.121
600	1.492	1.589	1.541	1.421
300	0.974	0.999	0.987	0.867
150	0.622	0.632	0.627	0.507
75	0.432	0.412	0.422	0.302
37.5	0.289	0.310	0.300	0.180
18.75	0.208	0.198	0.203	0.083
0	0.121	0.119	0.120	

LOD

4.68 pg/ml

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of human HLA-DRA in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:1	Average %	91
	Range %	84-98
1:2	Average %	91
	Range %	85-95
1:4	Average %	96
	Range %	92-102
1:8	Average %	103
	Range %	96-108

RECOVERY

The recovery of human HLA-DRA spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	92	88-96
EDTA plasma (n=4)	96	90-101