

PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%


Three samples of known concentration were tested twenty times on one plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

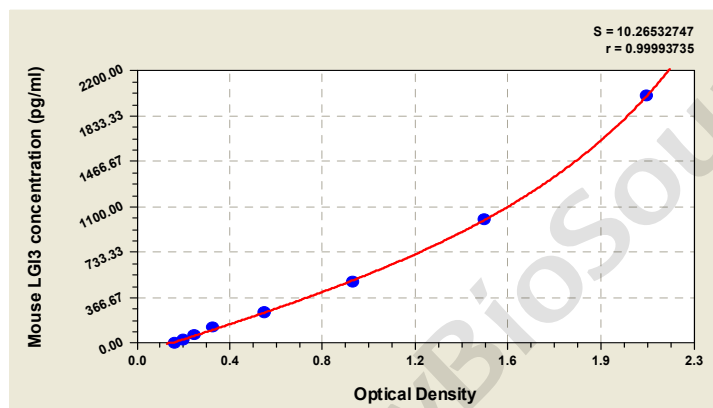

pg/ml	OD1	OD2	Average	Corrected	
2000	2.106	2.138	2.122	1.950	
1000	1.465	1.436	1.451	1.279	
500	0.895	0.918	0.907	0.735	
250	0.551	0.536	0.544	0.372	
125	0.329	0.331	0.330	0.158	
62.5	0.258	0.249	0.254	0.082	
31.25	0.201	0.210	0.206	0.034	
0	0.170	0.173	0.172		

LOD

7.8 pg/ml.

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of mouse LGI3 in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:1	Average %	95
	Range %	90-100
1:2	Average %	87
	Range %	83-90
1:4	Average %	102
	Range %	98-108
1:8	Average %	89
	Range %	85-94

RECOVERY

The recovery of mouse LGI3 spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	108	104-112
EDTA plasma (n=4)	100	95-105