

## **PRECISION**

**Intra-assay Precision (Precision within an assay): CV%<8%**


Three samples of known concentration were tested twenty times on one plate to assess.

**Inter-assay Precision (Precision between assays): CV%<10%**

Three samples of known concentration were tested in twenty assays to assess.

## **TYPICAL DATA**

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

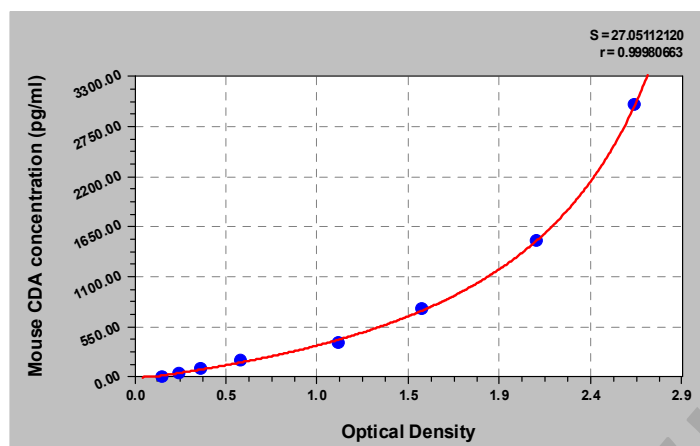

pg/ml	OD1	OD2	Average	Corrected
3000	2.711	2.578	2.645	2.485
1500	2.154	2.102	2.128	1.968
750	1.530	1.515	1.523	1.363
375	1.072	1.102	1.087	0.927
187.5	0.573	0.570	0.572	0.412
93.75	0.356	0.373	0.365	0.205
46.88	0.239	0.251	0.245	0.085
0	0.158	0.161	0.160	

## **LOD**

11.72pg/ml

**LINEARITY**

To assess the linearity of the assay, samples were spiked with high concentrations of mouse CDA in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:1	Average %	97
	Range %	93-101
1:2	Average %	85
	Range %	80-90
1:4	Average %	91
	Range %	85-95
1:8	Average %	94
	Range %	90-98

**RECOVERY**

The recovery of mouse CDA spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	106	101-113
EDTA plasma (n=4)	93	88-97