

PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%


Three samples of known concentration were tested twenty times on one plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

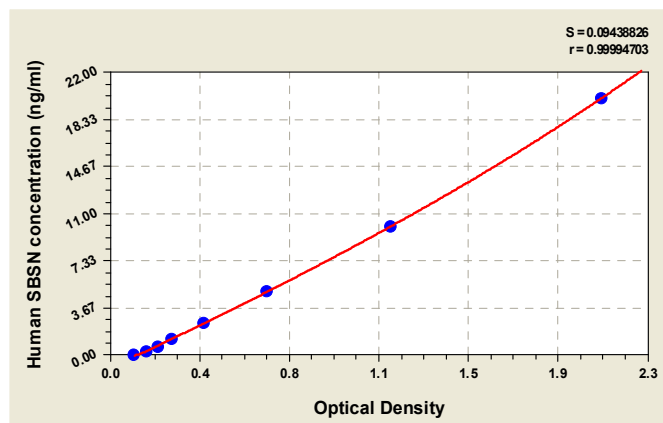

ng/ml	OD1	OD2	Average	Corrected
20	2.012	2.158	2.085	1.973
10	1.174	1.216	1.195	1.083
5	0.655	0.689	0.672	0.560
2.5	0.412	0.403	0.408	0.296
1.25	0.281	0.261	0.271	0.159
0.625	0.208	0.211	0.210	0.098
0.312	0.162	0.169	0.166	0.054
0	0.111	0.112	0.112	

LOD

0.078 ng/ml

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of human SBSN in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:1	Average %	97
	Range %	92-103
1:2	Average %	88
	Range %	83-94
1:4	Average %	92
	Range %	87-98
1:8	Average %	101
	Range %	93-104

RECOVERY

The recovery of human SBSN spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	87	81-93
EDTA plasma (n=4)	102	97-106