

PRECISION

Intra-assay Precision (Precision within an assay): CV%<8%


Three samples of known concentration were tested twenty times on one plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

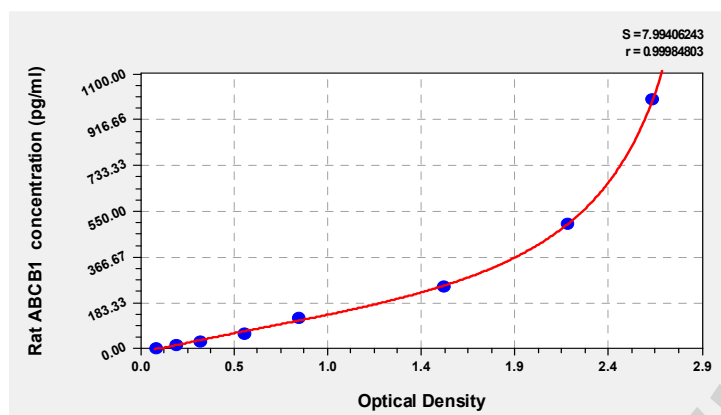

pg/ml	OD1	OD2	Average	Corrected
1000	2.655	2.546	2.601	2.511
500	2.141	2.200	2.171	2.081
250	1.556	1.533	1.545	1.455
125	0.777	0.845	0.811	0.721
62.5	0.532	0.545	0.539	0.449
31.2	0.325	0.302	0.314	0.224
15.6	0.188	0.188	0.188	0.098
0	0.089	0.090	0.090	

LOD

3.9 pg/ml

LINEARITY

To assess the linearity of the assay, samples were spiked with high concentrations of rat ABCB1 in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:1	Average %	88
	Range %	81-95
1:2	Average %	91
	Range %	86-97
1:4	Average %	97
	Range %	91-102
1:8	Average %	102
	Range %	97-107

RECOVERY

The recovery of rat ABCB1 spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	103	95-110
EDTA plasma (n=4)	93	85-98