

## **PRECISION**

**Intra-assay Precision (Precision within an assay): CV%<8%**


Three samples of known concentration were tested twenty times on one plate to assess.

**Inter-assay Precision (Precision between assays): CV%<10%**

Three samples of known concentration were tested in twenty assays to assess.

## **TYPICAL DATA**

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

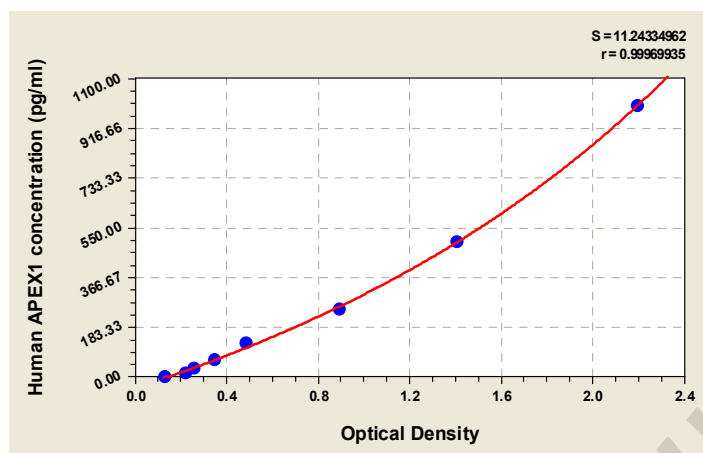

pg/ml	OD1	OD2	Average	Corrected
1000	2.253	2.155	2.204	2.058
500	1.398	1.434	1.416	1.270
250	0.916	0.900	0.908	0.762
125	0.507	0.487	0.497	0.351
62.5	0.351	0.366	0.359	0.213
31.25	0.268	0.272	0.270	0.124
15.6	0.229	0.234	0.232	0.086
0	0.147	0.145	0.146	

## **LOD**

3.9 pg/ml.

**LINEARITY**

To assess the linearity of the assay, samples were spiked with high concentrations of human APEX1 in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

	Sample	Serum(n=4)
1:1	Average %	104
	Range %	93-108
1:2	Average %	104
	Range %	95-107
1:4	Average %	91
	Range %	84-95
1:8	Average %	92
	Range %	84-95

**RECOVERY**

The recovery of human APEX1 spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

Sample Type	Average % Recovery	Range
Serum (n=5)	100	89-105
EDTA plasma (n=4)	104	93-108